


# HBsAg

One step qualitative sandwich immunochromatographic assay for detection of Hepatitis 'B' surface antigen

- Direct Sandwich Assay
- Specimen Type: Serum, Plasma
- Specimen Volume: 1 Drop (Dropper Provided) or 50  $\mu$ l
- Result Interpretation Time: 20 Minutes
- Storage Temperature: 2 to 30°C
- Analytical Sensitivity: 0.5 ng / ml
- Diagnostic Sensitivity: > 98.0 %
- Diagnostic Specificity: > 99.5%
- Long shelf life



## Test protocol

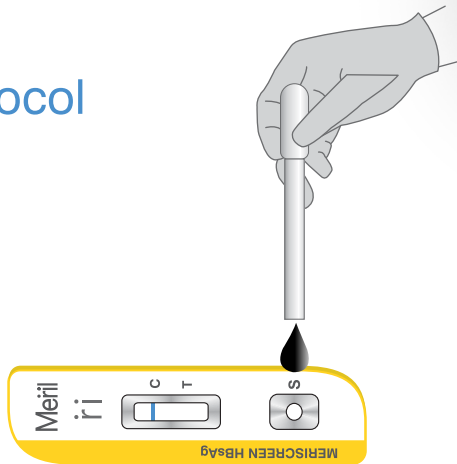


Add 1 drop or 50 $\mu$ l of plasma or serum into the sample well using a Dropper/micropipette



Read results at 20 minutes

## Interpretation of test results

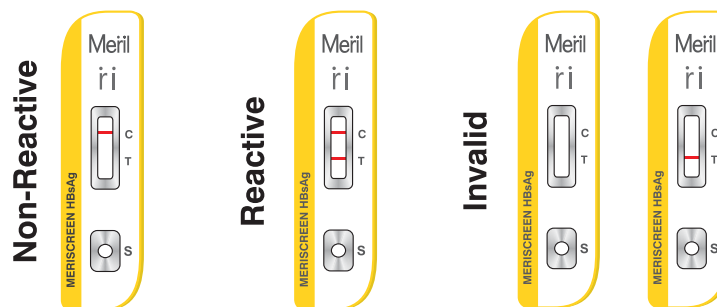

Product Code	Material Description	Pack Size
RPDHBV - 01	MERISCREEN HBsAg	50 Test